

# Artifact Removal in EEG Data

Maie LE GUILLY - Rasmus HAARSLEV - Troels THOMSEN

## 1 Purpose

Our purpose is to remove noise from EEG data, in order to allow for real-time signal analysis.

## 2 Methods

We use PCA to identify components with unusually large eigenvalues. If the eigenvalue is above the threshold found during calibration, there should be an artifact in that component and we try to remove it.

To compute the threshold for rejecting components we use a clean dataset as a training set. We calibrate the threshold on the first 20% of the dataset and then add simulated artifacts to the rest of it.

We use three different kind of threshold:

- The maximum eigenvalue found in the windows of the training set
- The average of the eigenvalues of the windows of the training set
- The average of the maximum of each window of the training set

## 3 Results

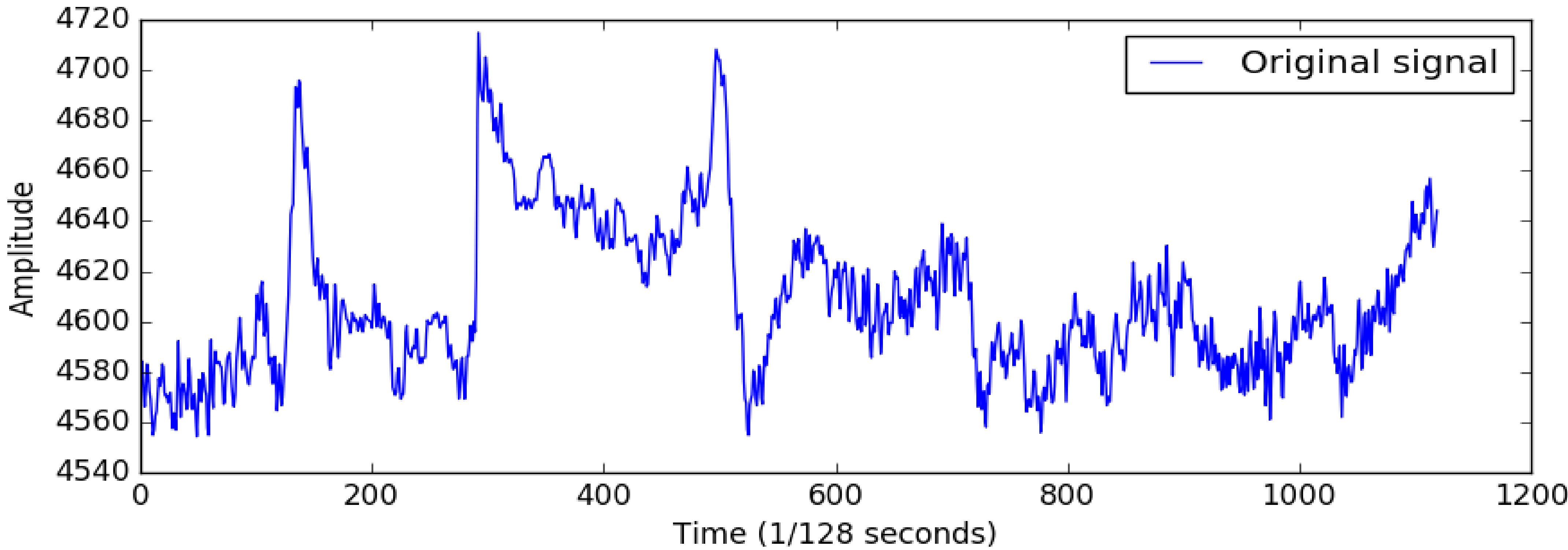


Figure 1: representation of one channel of an EEG dataset

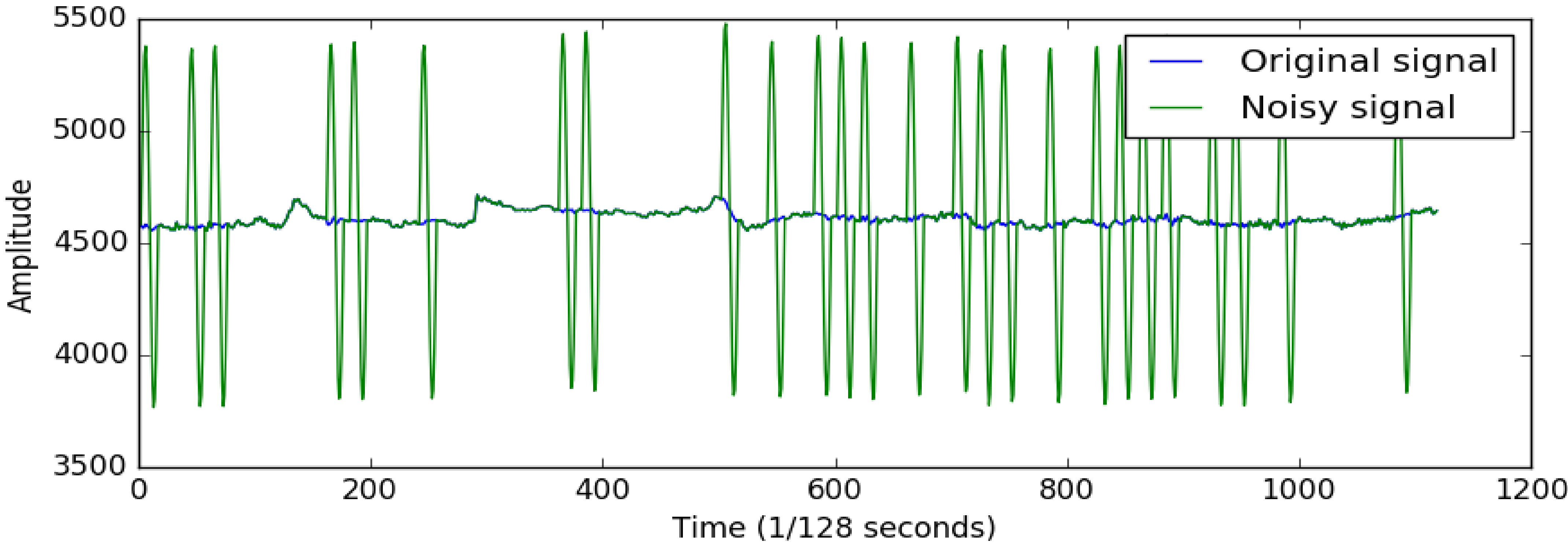


Figure 2: original dataset with simulated artifacts added

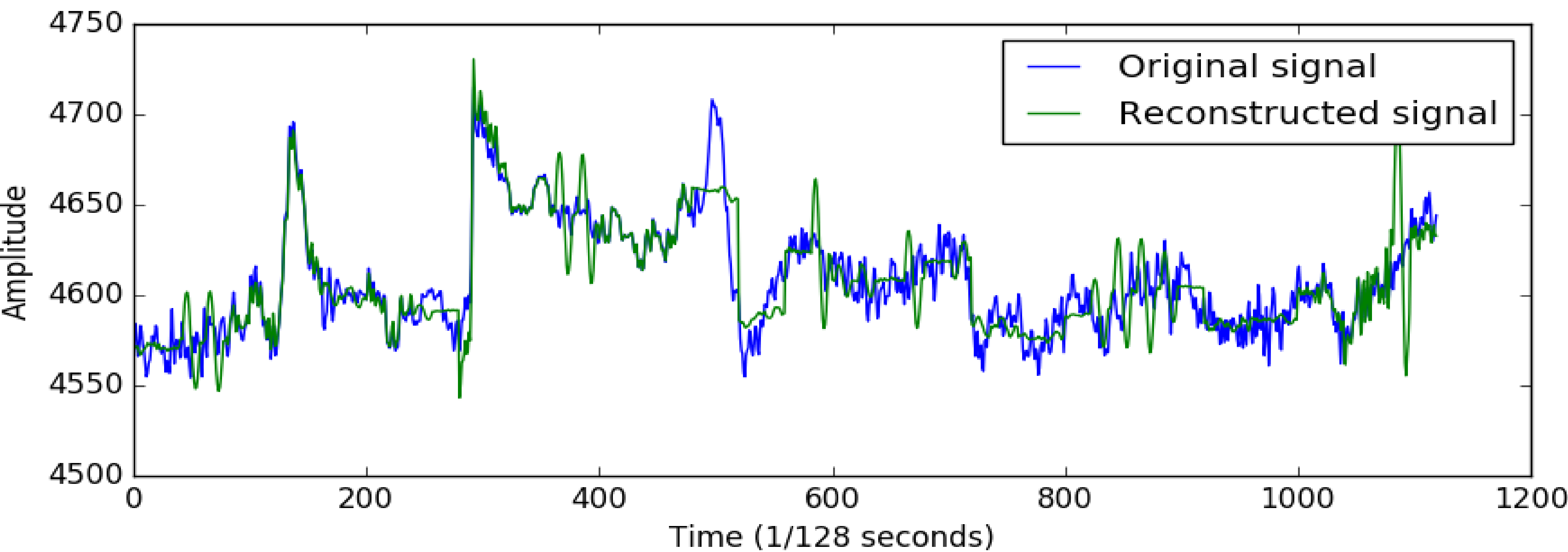


Figure 3: original dataset together with the dataset reconstructed from artifact dataset

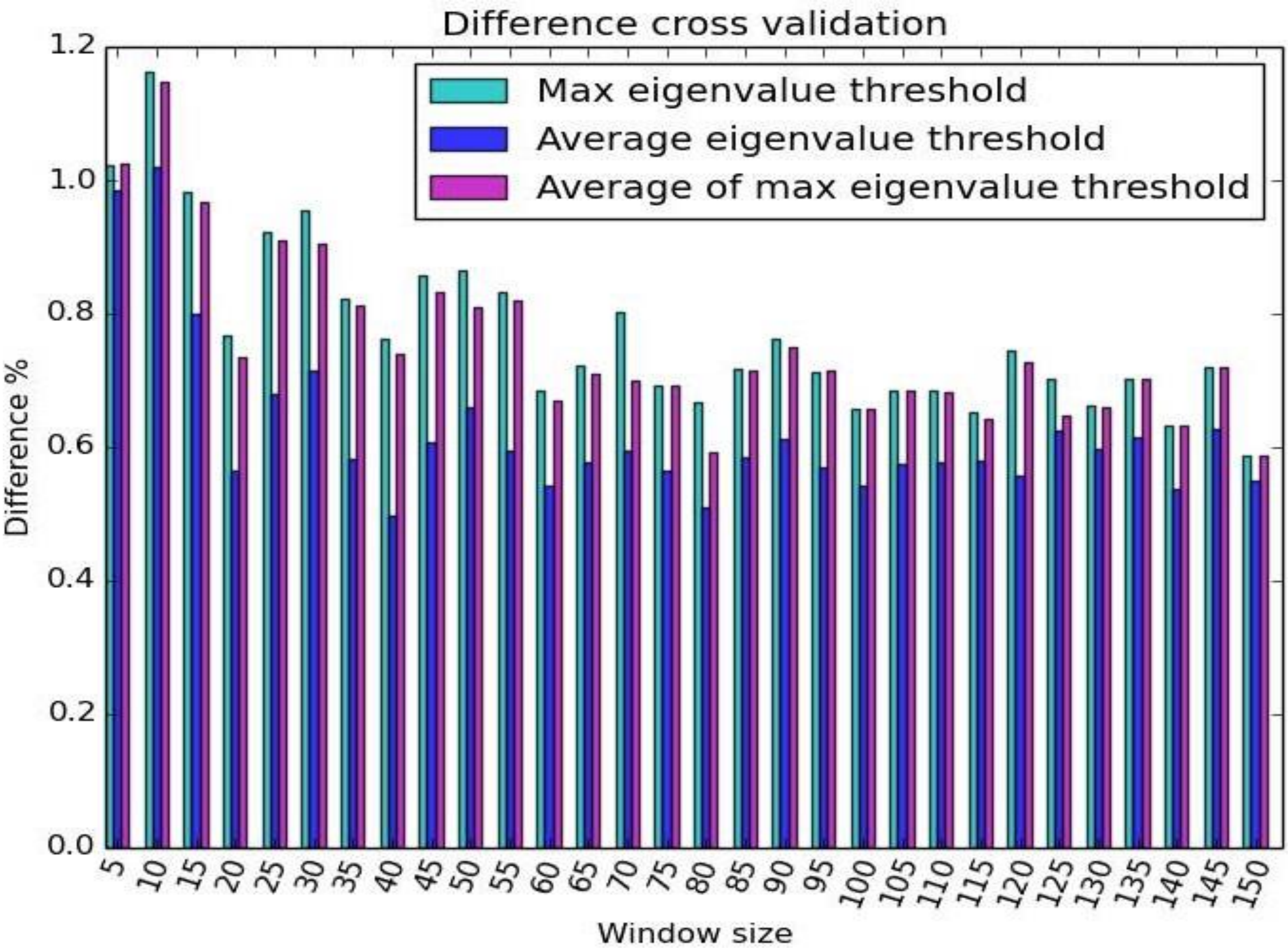


Figure 4: cross-validation between thresholds and window size on a 1400 sample dataset

Threshold	Max	AVG	MAX_AVG
Sensitivity	0.92	1.0	1.0
Specificity	0.22	0.0	0.16

Figure 5: Sensitivity and specificity for the different threshold

For the average threshold and a window size of 40 samples with randomly added artifacts we have an average loss in the reconstructed dataset of 0.52%. If no artifacts are added this becomes a 0.027% loss. For the same window size we can remove artifacts at a rate of 81Hz allowing real time artifact removal.

## 4 Conclusion

We can remove artifacts with an average loss of 0.52% at a rate estimated to support real time analysis.